


Reduction made easy

The JET-X Bar Ankle Spanning Kit - Multiple Pin Clamp offers a simple solution for the treatment of tibial fractures.

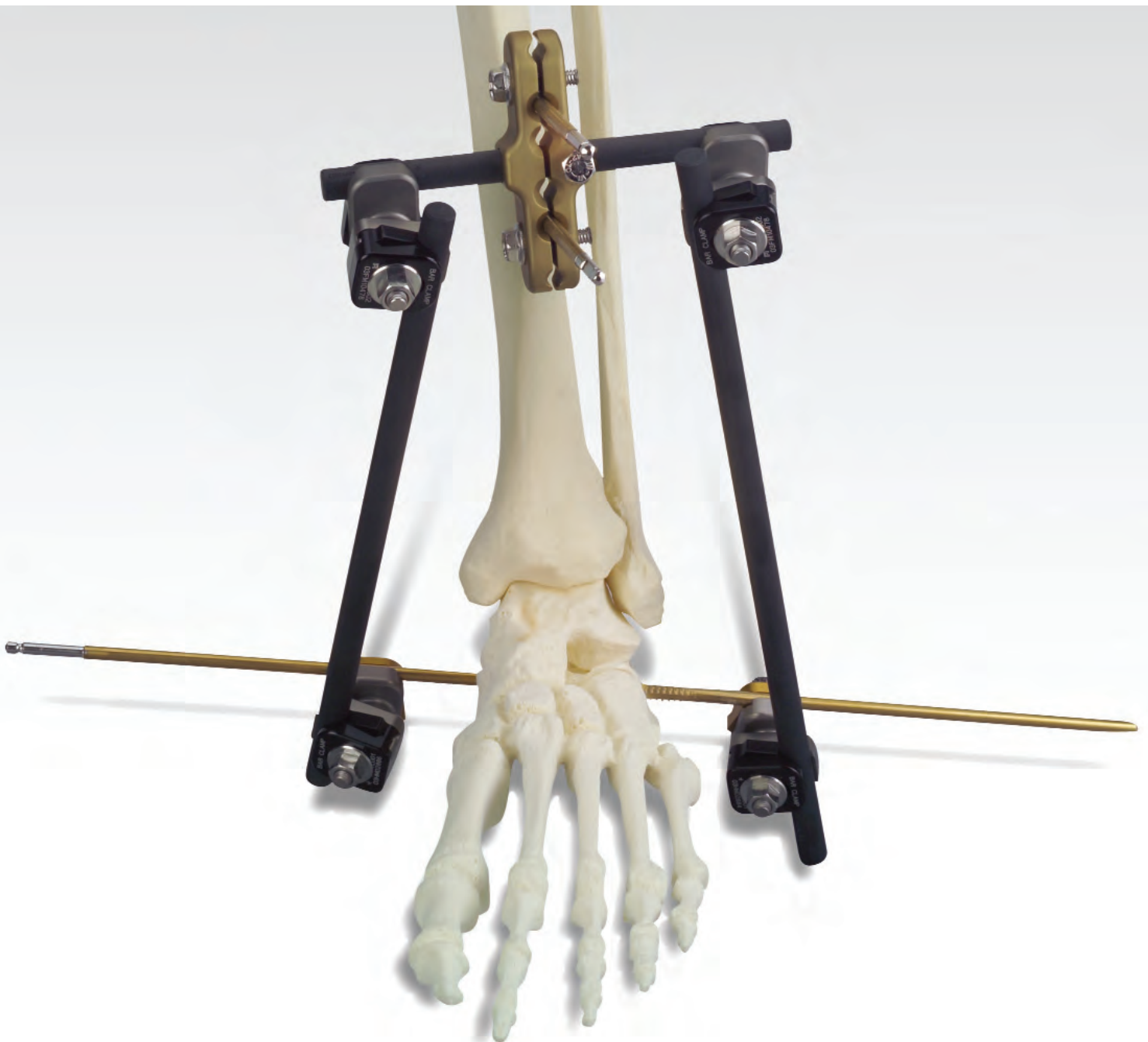

JET-X[◇] Bar Ankle Spanning Kit - Multiple Pin Clamp

The JET-X Bar Ankle Spanning Kit - Multiple Pin Clamp is designed for ease of application. The JET-X Bar to Pin Clamp and JET-X Bar to Bar Clamp facilitates fracture reduction while providing stability in a low profile and lightweight design. The ball joint feature provides the ability to easily position the ankle in multiple planes to aid in fracture reduction.

Design Features



Lightweight, radiolucent carbon fiber bars.

Ball joints angulate 25° allowing an extra degree of freedom to reduce fractures in multiple planes.

Clamps have independent locking mechanisms, allowing reduction after the frame is attached to the pins.

JET-X Bar Ankle Spanning Kit Multiple Pin Clamp

Catalog No. 7106-2701 includes:

Description	Qty.
Bar to Pin Clamp MR Safe	2
Bar to Bar Clamp MR Safe	2
Multiple Pin Clamp	1
100mm Bar	1
150mm Bar	2
5mm x 50mm TiN Traction Pin	2
5mm x 40mm TiN Half Pin	3
Protective Cap, Black – 5mm	4
Drill Sleeve with Trocar	1
3.5mm Step Drill/AP Conn.	1
Pin Driver with 10mm Hex	1
10mm Disposable Wrench	1



Sterile Kit for JET-X Bar Ankle Spanning Kit - Multiple Pin Clamp