

Construct guide



Table of contents

Ankle spanning (version 1)	2
Ankle spanning (version 2)	3
Traveling traction.....	4
Tibia.....	5
Femur	6
Knee spanning.....	7
Pelvis	8
Elbow	9
Distal radius	10
Trays.....	11
Sterile kits.....	12


JET-X^o Ankle spanning (version 1)

Description	Qty
10.5mm Bar	2
10.5mm to 5mm Pin Clamp	4
5mm Half Pin	2
5mm Traction Pin	1
10.5mm to 4mm Pin Clamp	1
4mm Half Pin	1



JET-X^o Ankle spanning (version 2)

Description	Qty
10.5mm Bar	2
6 Hole Multi Pin Clamp	1
Freedom Post	2
5mm Half Pin	2
10.5mm Bar to 5mm Pin Clamp	2
5mm Traction Pin	1
10.5mm to 4mm Pin Clamp	1
4mm Half Pin	1

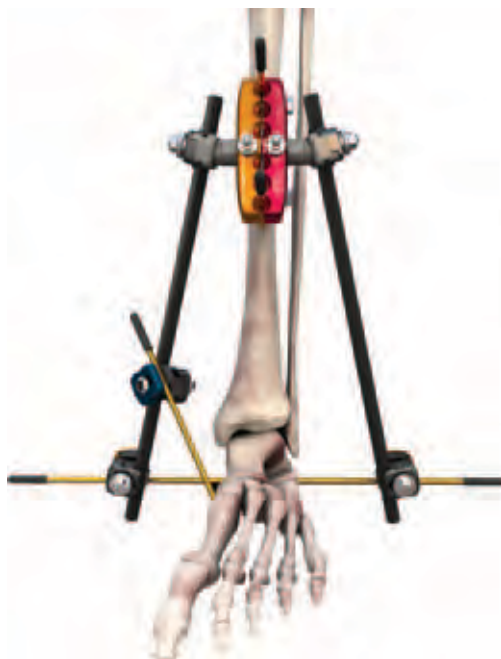

JET-X[®] Traveling traction

Description	Qty
10.5mm Bar	2
10.5mm to 5mm Pin Clamp (Freedom or Quick)	4
5mm Traction Pin	2




JET-X[®] Basic tibial shaft

Description	Qty
10.5mm Bar	1
10.5mm Bar to 5mm Pin Clamp	4
5mm Half Pin	4



JET-X^o Basic femur

Description	Qty
10.5mm Bar	1
10.5mm to Pin Clamp	4
Half Pins	4

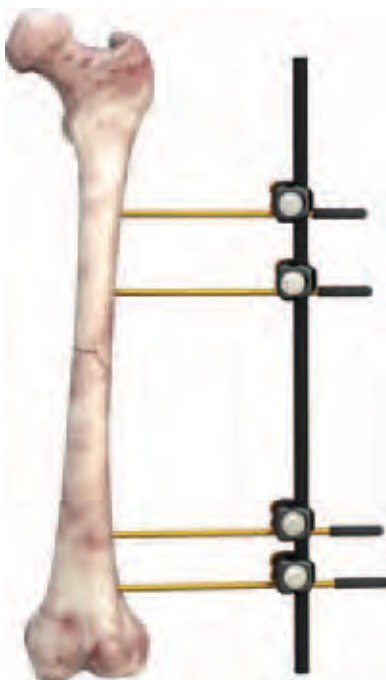

JET-X^o Knee spanning

Description	Qty
10.5mm Bar	3
10.5mm Bar to 10.5mm Bar Clamp	2
10.5mm Bar to Pin Clamp	4
Half Pins	4



JET-X^o Pelvic

Description	Qty
10.5mm Bar	2
10.5mm Bar to Bar Clamp	1
10.5mm Bar to Pin Clamp	2
Half Pins	2



JET-X^o Elbow spanning

Description	Qty
10.5mm Bar	3
10.5mm to 10.5mm Bar Clamp	3
10.5mm to 5mm Pin Clamp	2
10.5mm Bar to 4mm Pin Clamp	2
5mm Half Pin	2
4mm Half Pin	2



JET-X[®] Distal radius

Description	Qty
6mm Bar	1
Mini Double Pin Clamp	2
4mm Half Pin	4



JET-X[®] Trays

**Large Clamp
Module Tray**



**Large Partial Half Pin
and Instrument
Module Tray**



Large Bar Module Tray



**JET-X Mini Instruments
and Trays**



JET-X[®] Kits

**Ankle Spanning Kit,
Standard Central Body
with TiN Half Pins,
Sterile**

Cat. No. 7105-1701



**Ankle Spanning Kit,
Multiple Pin Clamp with
TiN Half Pin, Sterile**

Cat. No. 7106-2701



**Ankle Spanning Kit,
Single Pin Clamp with
TiN Half Pin, Sterile**

Cat. No. 7106-2702



**Travel Traction Kit,
600mm Bar with TiN
Traction Pin, Sterile***

Cat. No. 7106-2703



**Tibia Kit, 350mm Bar
with TiN Half Pin,
Sterile***

Cat. No. 7106-2705



**JET-X Mini Bar External
Fixator Kit, Sterile**

Cat. No. 7106-4701



Orthopaedics

Smith & Nephew, Inc.
1450 Brooks Road,
Memphis, TN 38116, USA

www.smith-nephew.com

Telephone: 1-901-396-2121
Information: 1-800-821-5700,
Orders and Inquiries: 1-800-238-7538

°Trademark of Smith & Nephew. Reg. US Pat. & TM Off.
All Trademarks acknowledged. ©2009 Smith & Nephew, Inc.

7108-0877 REV0 10/09